


Advanced technology for  
transforaminal injections

**NeuroTherm<sup>®</sup>**  
**TRUCATH<sup>◊</sup>**  
Spinal Injection System



Cervical transforaminal injections carry an inherent risk of serious complications. The TRUCATH® Spinal Injection System is an innovative device specifically designed to keep the needle tip positioned away from sensitive structures.

- Blunt, flexible catheter designed to deflect off of vasculature and neural elements
- Side injection orifice for injection of contrast and therapeutic agents proximal of the catheter tip
- Visual and tactile direction indicator for desired orientation of needle bevel and deployment of catheter
- Radiopaque tantalum marker within the catheter tip for fluoroscopic visualization
- Microbore tubing with low flush volume

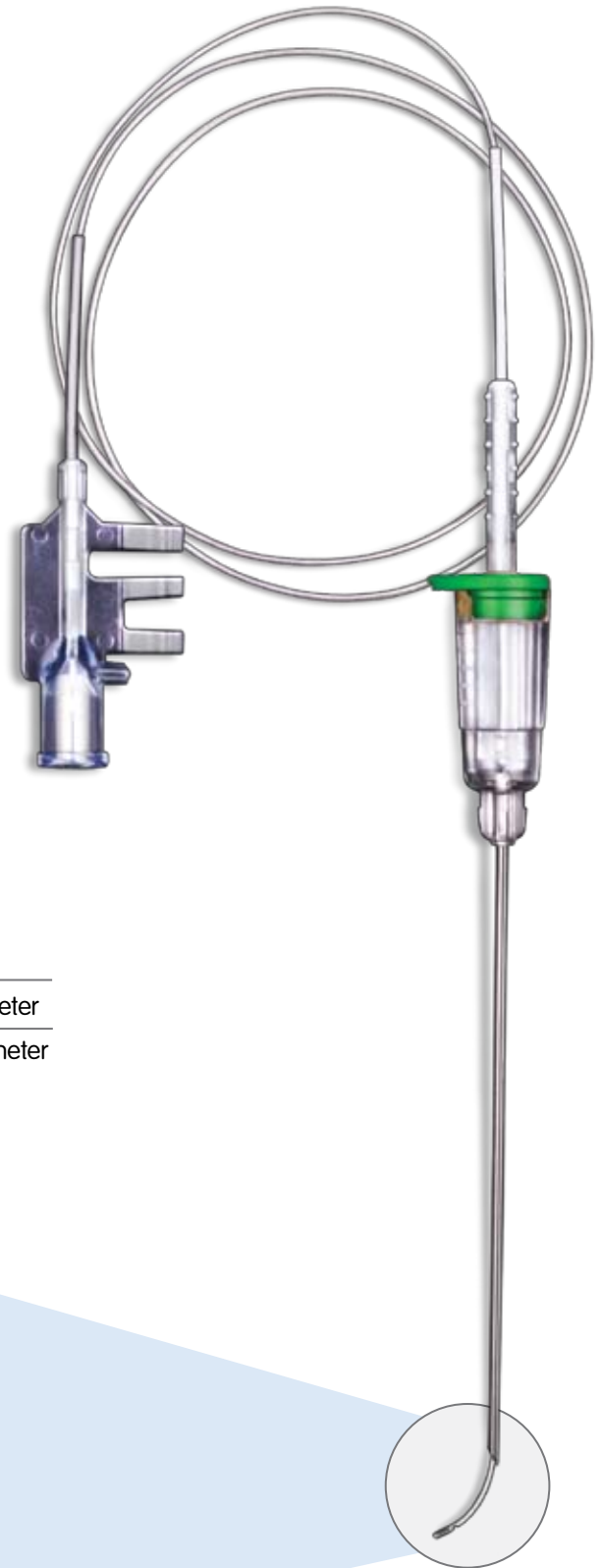

### Ordering information:

Item number	Description
NT7210320	TRUCATH Spinal Injection System, 22G x 1.8"(46mm), 25G Catheter
NT7210321	TRUCATH Spinal Injection System, 20G x 4.5"(114mm), 23G Catheter

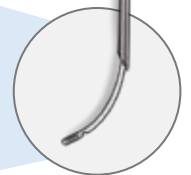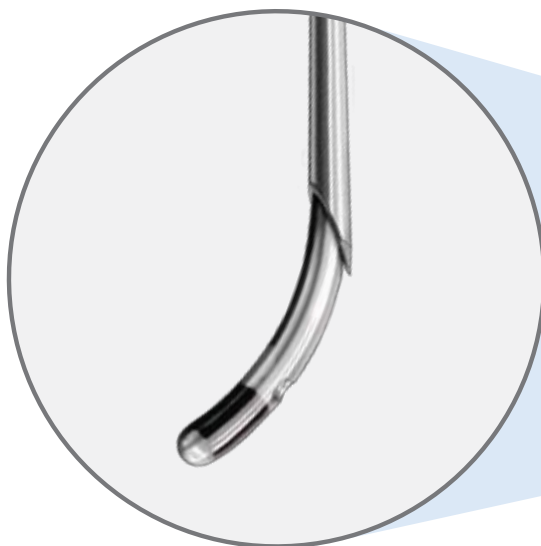

**Distributed by:**  
**NeuroTherm®**  
 30 Upton Drive, Suite 2  
 Wilmington, MA 01887, United States  
 1-888-655-3500  
 1-978-658-2378 FAX  
 usasales@neurotherm.com